



## CONSUMER PRODUCTS

### TABLE OF CONTENTS

- ◆ GUARDIAN Forehead Thermometers, Bulk
- ◆ FEVERSCAN ULTRA Forehead Thermometers, Packaged
- ◆ EAR PUTTY Silicone Ear Plugs
- ◆ CLEAN FIT™ Foam Ear Plugs
- ◆ PROVITA Smoke Hood

**Tecfen Corporation**

5860-C Hollister Avenue  
Santa Barbara, CA 93117 U.S.A

Tel: (805) 967 1153  
Fax: (805) 967 1295

Home: [www.tecfen.com](http://www.tecfen.com)  
Email: [tecfen@tecfen.com](mailto:tecfen@tecfen.com)

# Guardian Forehead Thermometer, BULK



Taking your temperature has never been so clean and easy

- Does not break
- Easy to use
- Fast - Indicates in 15 seconds
- Safe - No risk of Mercury spillage
- Easy to read for young and old
- Easy to carry



The "DIGITAL" Incorporates the benefits of previous forehead thermometer technology together with an extra advantage; a more precise display of body temperature is indicated in 1/2 degree temperature increments. Within 15 seconds the core temperature of the individual is indicated in one of nine easy-to-read windows. Other attractive features includes its comfortable, unbreakable design and convenient storage capabilities.

**THE DIRECTIONS ARE EASY:** Make sure the forehead is clean and dry. Firmly apply the thermometer onto the middle of the forehead, holding it at both ends, without touching the numbers. Wait until colors stop changing, usually 15-20 seconds. For the correct temperature, READ THE GREEN COLOR ONLY. Disregard Blue and Tan colors. Do not use the thermometer directly exposed to sunlight or very bright lamps.

**THE TECHNOLOGY:** The Digital Forehead Thermometer is a quick, convenient, and accurate thermometer using Liquid Crystal Technology. Liquid Crystals are the most thermally sensitive materials known, possessing a unique state of matter: the optical properties of both liquids and crystals. The colorful diffracting liquid crystals are cholesteric esters which change colors when their molecules, which are in a spiral structural form, are caused to change due to temperature variations. Such molecular movement results in visible color shifts which are evident to the eye.

Research has shown that forehead temperatures closely follow internal, or "core" temperature movements. The central forehead is notably free of fat or obscuring and insulating tissue masses. Digital Thermometers when placed on the forehead correspond directly to core temperature variations. Measuring from the surface allows you to quickly and easily determine any temperature trends.



**Tecfen Corporation**

5860-C Hollister Avenue  
Santa Barbara, CA 93117 U.S.A




Tel: (805) 967 1153  
Fax: (805) 967 1295



**Home:** [www.tecfe.com](http://www.tecfe.com)  
Email: [tecfe@tecfe.com](mailto:tecfe@tecfe.com)

# *FEVERSCAN ULTRA*

## *Forehead Thermometer*

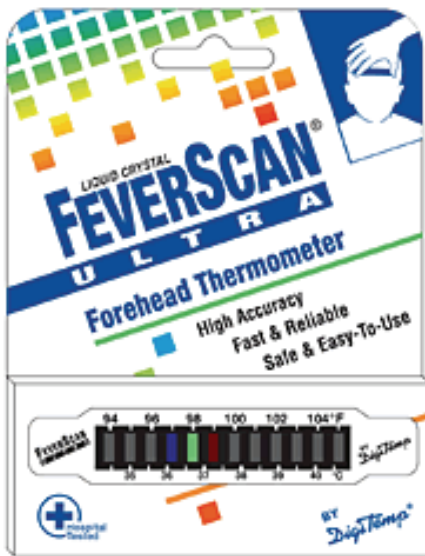

**Forehead Thermometer Monitor**

**Reliable replace back of hand**

**Alerts to possible fever**

**Easy to Read**

**Fast, Latex Free**


**Temperature Range:**

**94-105°F (34.5-40.5°C)**

**Accurate to  $\pm 1^{\circ}\text{F}$  ( $\pm 0.6^{\circ}\text{C}$ )**

**Mercury Free-No Glass**

**Reusable**

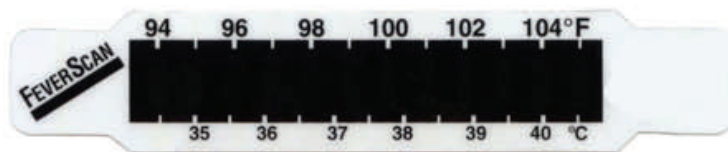


5860-C Hollister Avenue  
Santa Barbara, CA 93117 U.S.A.

Telephone: 805 967-1153  
FAX: 805 967-1295

Website: [www.tecfe.com](http://www.tecfe.com)  
email: [tecfe@tecfe.com](mailto:tecfe@tecfe.com)

# EAR PUTTY

## Silicone Ear Plugs

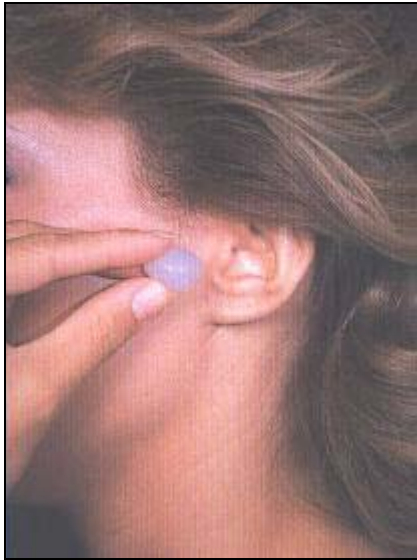

Use for:

- Swimming • Bathing • Showering • Sleeping • Studying
- Traveling

Made from silicone, a polymer that is hypoallergenic, nontoxic, and water-repellent.

Noise reduction rating of 20 decibels.

- Effective
- Easy to use
- Reusable
- Soft Silicone
- Comfortable

WATERPROOF - BLOCKS NOISE



SB-30024



Packaged in shrink-wrap  
24 cases per display box



2 plugs per small  
carrying case



5860-C Hollister Avenue  
Santa Barbara, CA 93117 U.S.A

Tel: (805) 967 1153  
Fax: (805) 967 1295

Home: [www.tecfe.com](http://www.tecfe.com)  
Email: [tecfe@tecfe.com](mailto:tecfe@tecfe.com)

**SB-30108**

Packaged 8 plugs  
with a small plastic carrying case,  
on a display card



## WATERPROOF - BLOCKS NOISE

Note: Discard Ear Putty when dirty or no longer sticky.



5860-C Hollister Avenue  
Santa Barbara, CA 93117 U.S.A



Tel: (805) 967 1153  
Fax: (805) 967 1295



Home: [www.tecfe.com](http://www.tecfe.com)  
Email: [tecfe@tecfe.com](mailto:tecfe@tecfe.com)

# CLEAN FIT™ Foam Ear Plugs



**30400 Dirt-Resistant Surface  
Reduced Annoying Sounds**

Smooth, dirt resistant, non-allergenic surface. Ideal for working and operating loud machinery.

Noise Reduction Rating (NRR): 29.

- Safe
- Easy to Use
- Effective
- Self-Adjusting
- Reusable
- Soft
- Comfortable



Packaged 20 plugs on a display card.



5860-C Hollister Avenue  
Santa Barbara, CA 93117 U.S.A.

Telephone: 805 967-1153  
FAX: 805 967-1295

Website: [www.tecfe.com](http://www.tecfe.com)  
email: [tecfe@tecfe.com](mailto:tecfe@tecfe.com)

# *PROVITA Smoke Hood*



## **Life-Saving Smoke Hood Provides Vital Minutes**

More than 80% of fire-related deaths actually are caused by smoke and chemically-poisoned air. Give yourself and family 20 minutes of evacuation time with this smoke hood - its transparent, heat-resistant shield is gas-impermeable, with a 360 degree charcoal filter and a rubber seal at the neck.

- Lightweight and Compact  
*Weight: 160 grams (6 oz.)*  
*Package Dimensions: 20cm x 15cm x 2cm (8" x 6" x 1")*
- Sealed package can easily be kept in briefcase, desk drawer, or handbag
- One size fits all
- Proven to protect from smoke and 98% of toxins created by fire, including  
*Chlorine*  
*Ammonia*  
*Hydrogen chloride*  
*Hydrogen sulfide*  
*Acrolein*



### **Tecfen Corporation**

5860-C Hollister Avenue  
Santa Barbara, CA 93117 U.S.A



Tel: (805) 967 1153  
Fax: (805) 967 1295

Home: [www.tecfe.com](http://www.tecfe.com)  
Email: [tecfe@tecfe.com](mailto:tecfe@tecfe.com)