



MEDICAL SUNDRIES

TABLE OF CONTENTS

- ◆ Digital Forehead Thermometers
- ◆ Tecfen All Weather Blanket
- ◆ Tecfen Disposable Blanket

Digital Forehead Thermometer



Taking your temperature has never been so clean and easy

- Does not break
- Easy to use
- Fast - Indicates in 15 seconds
- Safe - No risk of Mercury spillage
- Easy to read for young and old
- Easy to carry



The "DIGITAL" Incorporates the benefits of previous forehead thermometer technology together with an extra advantage; a more precise display of body temperature is indicated in 1/2 degree temperature increments. Within 15 seconds the core temperature of the individual is indicated in one of nine easy-to-read windows. Other attractive features includes its comfortable, unbreakable design and convenient storage capabilities.

THE DIRECTIONS ARE EASY: Make sure the forehead is clean and dry. Firmly apply the thermometer onto the middle of the forehead, holding it at both ends, without touching the numbers. Wait until colors stop changing, usually 15-20 seconds. For the correct temperature, READ THE GREEN COLOR ONLY. Disregard Blue and Tan colors. Do not use the thermometer directly exposed to sunlight or very bright lamps.

THE TECHNOLOGY: The Digital Forehead Thermometer is a quick, convenient, and accurate thermometer using Liquid Crystal Technology. Liquid Crystals are the most thermally sensitive materials known, possessing a unique state of matter: the optical properties of both liquids and crystals. The colorful diffracting liquid crystals are cholesteric esters which change colors when their molecules, which are in a spiral structural form, are caused to change due to temperature variations. Such molecular movement results in visible color shifts which are evident to the eye.

Research has shown that forehead temperatures closely follow internal, or "core" temperature movements. The central forehead is notably free of fat or obscuring and insulating tissue masses. Digital Thermometers when placed on the forehead correspond directly to core temperature variations. Measuring from the surface allows you to quickly and easily determine any temperature trends.



Tecfen Corporation

5860-C Hollister Avenue
Santa Barbara, CA 93117 U.S.A



Tel: (805) 967 1153
Fax: (805) 967 1295



Home: www.tecfe.com
Email: tecfe@tecfe.com

Tecfen All Weather Blanket




"The all reason – all season outdoor blanket"

This versatile Outdoor Safety Blanket provides ultimate personal protection from the harsh elements of nature. Its unique combination of metal and specific plastic films form a perfect barrier against wind and moisture penetration. The blanket radiates the transmission of body heat to keep the individual warm and dry. Used by militaries, police, disaster relief agencies, fire & medical organizations worldwide.

Uses:

- Ground Cover
- Emergency Shelter
- Ground to Air Signal
- Heat Reflector
- Medical Blanket
- Boating Blanket
- Sleeping Bag Liner
- Improvised Cooler
- Desert Sunshade
- Picnic Blanket
- Stadium Blanket

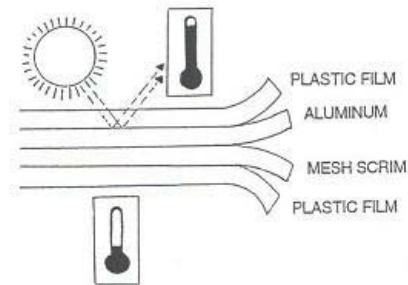

Specifications:

Color: Red/Silver

Size: 152 cm x 213 cm

Weight: 340 grams

Package Size: 29 cm x 20 cm x 2 cm



Tecfen Corporation

5860-C Hollister Avenue
Santa Barbara, CA 93117 U.S.A

Tel: (805) 967 1153
Fax: (805) 967 1295



Home: www.tecfen.com
Email: tecfen@tecfen.com

Tecfen Disposable Blanket

TE-105 TECFEN DISPOSABLE BLANKET



The Disposable Blanket is highly absorbent, soft and pliable blanket made of 2-ply bleached non-wet strength paper tissue laminated onto a thick yellow low density polyethylene film. The blanket provides excellent protection, padding, and cushioning. In emergency situations, it is used because of its high absorption of body fluids.

Uses:

- Non-sterile emergency medical use
- Protective underlay for many different uses
- Emergency Blanket for body heat retention
- Assists in prevention of post-accident shock
- Emergency Tarp for weather protection
- Reversible for liquid absorption or water repellence
- Soft pliable material
- Light weight & compact packaging

Specifications:

Color: Yellow
Size: 137cm x 203cm
Weight: 180 grams



Tecfen Corporation

5860-C Hollister Avenue
Santa Barbara, CA 93117 U.S.A

Tel: (805) 967 1153
Fax: (805) 967 1295



Home: www.tecfe.com
Email: tecfe@tecfe.com